

## In-Sight™ Multi-Band Ligator

### Instructions for Use

 Consult Instructions for Use  Do Not Reuse

 Non-Sterile **Rx only**

 Do Not Use If Package Is Damaged

### IMPORTANT INFORMATION

Caution: Federal law restricts this device to sale by or on the order of a physician. Read all instructions carefully before use. They contain essential information on using this device safely and effectively. Keep these instructions in a safe, accessible location, as you may need to refer to them again. If you have any questions or comments about any information in these instructions, please contact Micro-Tech USA.

### INTENDED USE

The In-Sight Multi-Band Ligator is used to ligate internal hemorrhoids in adult patients only.

### CONTRAINDICATIONS

- Contraindications specific to hemorrhoidal banding include but are not limited to; severe inflammatory bowel disease, coagulopathy or anal strictures.

### POTENTIAL COMPLICATIONS

- Potential complications associated hemorrhoidal banding include but are not limited to, pain, bleeding, urinary symptoms, swelling, edema, tissue ulceration, and band dislodgement.
- Hemorrhoidal banding may result in severe pain if the procedure is performed below the dentate line.

### WARNINGS AND PRECAUTIONS

- The product is intended for single use only! DO NOT re-use, and/or reprocess. Re-use, or reprocessing may compromise the structural integrity of the device and/or lead to device failure which, in turn, may result in patient injury, illness or death. Re-use, or reprocessing may also create a risk of contamination of the device and/or cause patient infectious disease(s). Contamination of the device may lead to injury, illness or death of the patient. Micro-Tech assumes no liability with respect to instruments reused, or reprocessed.
- Do not use this instrument for any purpose other than its intended use.
- The product is only intended for adult populations.
- Hemorrhoid banding may result in severe pain if the procedure is performed below the dentate line.
- Banding should begin at the most proximal location from the anal sphincter and proceed distally because passing the ligator over a previously

placed band may dislodge band.

- The instrument is intended for use under the direct supervision of a suitably trained physician only. A thorough understanding of the technical principles, clinical applications, and associated risks is expected before usage.
- After use, this product may be a potential biohazard. Handle and dispose of in accordance with hospital, local and administrative laws and regulations.

**【 Product Name 】** In-Sight™ Multi-Band Ligator


**【 Production Date 】** See packaging

**【 Shelf Life 】** 1 year

This device is not made with natural rubber latex.

### STRUCTURE

The In-Sight Multi-Band Ligator provides illumination, suction, direct visualization and ligation capabilities in a single use device.

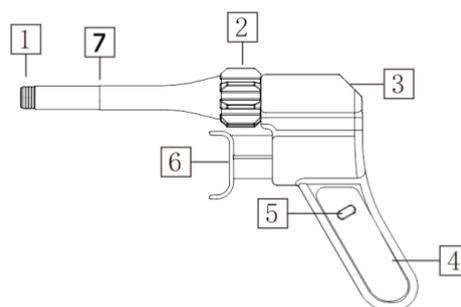

**Fig. 1 Diagram of In-Sight Multi-Band Ligator**

- Ligator bands
- Band actuation knob
- Viewing pathway
- Handle
- Light on/off button
- Vacuum drawing trigger
- Depth indicator

### PREPARATION

- Inspect the package before use for any damage. Do not use if package is damaged.
- Verify the expiration date. Do not use if expired.
- Open the package and carefully remove the ligator. Inspect its integrity and proper placement of bands. Confirm function of the light by pressing the light on/off button.
- If this device shows any signs of damage, or the light fails to illuminate, do not use. Do not attempt to repair a nonfunctional or damaged device.

### INSTRUCTIONS FOR USE

- Lubricate the anus, perform rectal exam and insert an anoscope.
- Activate the light of the ligator by pressing the on/off button on the ligator handle.
- Inspect the lumen of the distal rectum for hemorrhoids with the ligator.
- If ligation is indicated, grasp the ligator. Place the index finger on the vacuum drawing trigger and three other fingers and thumb around the ligator handle.
- Squeeze the vacuum drawing trigger to force air out of the device.
- Insert the ligator through the anoscope toward the hemorrhoids.

7. When the view through the anoscope becomes obstructed by the ligator, use the viewing pathway of the ligator to maintain direct visualization of the hemorrhoids.
8. Position the tip of the ligator just above the hemorrhoids and press gently against the targeted area.
9. Release pressure on the trigger to activate suction.
10. Inspect the amount of tissue sucked into the tip of the ligator.
11. Confirm that the patient feels pressure and not pain.
12. Place the tip of the index finger on the band actuation knob and rotate clockwise. The rotation of the actuation knob will deploy a single band off the distal tip of the ligator.
13. Remove the ligator from the anoscope for repeat inspection of banding or ligate additional tissue without pulling the anoscope out (repeat above steps as needed).
14. If direct visualization cannot be achieved, use the indicator groove on the ligator barrel to confirm tip location.
15. After the procedure is complete, turn off the light and dispose of the device in an approved manner.

**Note:** Repeat rectal examination should be performed after the final banding to document that the band was properly placed and that only mucosa was banded. Banded tissue must have a narrow "neck."

**Note:** If the patient experiences pain or if the banded tissue has a thick "neck" then pushing the band up from the base of the banded tissue must be performed to release deeper layers of the rectal wall.

### **STORAGE**

The product should be stored in a cool, dry, clean, well-ventilated, non-corrosive gas environment. Do not expose the package to organic solvent, ionizing radiation or ultraviolet radiation.

### **PRODUCT DISPOSAL**

After use, this product may be a potential biohazard. Handle and dispose of in accordance with hospital, local and administrative laws and regulations.

### **LIMITED WARRANTY AND DISCLAIMERS**

- Limited Warranty to Buyer. Micro-Tech USA warrants to Buyer that, for the earlier of one (1) year from the date of purchase, or until the product is used by Buyer, the products will be free from defects in materials and workmanship when stored and used in accordance with the instructions for storage and use provided by Micro-Tech USA and in accordance with applicable regulatory requirements. Descriptions or specifications appearing in Micro-Tech USA's literature are meant to generally describe the products and do not constitute any express warranties. In the event that Micro-Tech USA gives technical advice with respect to the product, it is agreed that such advice is given without any liability on Micro-Tech USA's part.

Any guarantee of specific properties of or in the products shall only be effective if and to the extent specifically confirmed by Micro-Tech USA in writing. These warranties shall not apply for product failure or deficiency due to improper storage, alteration, or the consequences of uses for which the products were not designed or that adversely affect the products' integrity, reliability, or performance.

- Disclaimer and Release. The warranties, obligations, and liabilities of Micro-Tech USA as set forth herein are exclusive. Buyer hereby waives all other warranties, express or implied, arising by law or otherwise, with respect to the products and any other goods or services delivered by buyer, including, but not limited to: (1) all other express and implied warranties, including any implied warranty of merchantability or fitness for a particular purpose, and (2) any implied warranty arising from course of performance, course of dealing, or usage of trade.
- Implied Warranties. The purchase of products may be subject to laws in the territories applicable to the sale of the products by Micro-Tech USA to Buyer that impose implied warranties, conditions, or obligations that cannot be excluded, restricted, or modified, or can only be excluded, restricted, or modified to a limited extent. The provisions of Paragraphs 2 and 4 shall apply to the greatest extent allowed by such laws.
- Limitation of Liability. Except to the extent prohibited by applicable law, Micro-Tech USA's liability under this warranty is limited to: (a) the replacement of the products or the re-supply of equivalent products; (b) the repair of the products or payment of the cost of repairing the products; or (c) payment of the cost of replacing the products or acquiring equivalent products. Micro-Tech USA shall have no obligation or liability, whether arising in contract (including warranty), tort (including active, passive, or imputed negligence, strict liability, or product liability) or otherwise, for any special, consequential, punitive, incidental, or indirect damages, or for loss of use, loss of revenue, loss of business, lost profit, or other financial loss arising out of or in connection with any product or other goods or services furnished by Micro-Tech USA, even if Micro-Tech USA was aware of the possibility of such damages or loss.



**Manufactured for  
Micro-Tech USA Inc.**  
2855 Boardwalk Drive, Ann Arbor  
MI 48104 USA  
Tel: 734-259-3768  
Toll free: 877-552-4027  
[info@micro-tech-usa.com](mailto:info@micro-tech-usa.com)  
Made in P.R. China

In-Sight™ is a trademark of Micro-Tech Endoscopy.

MDHL100000 Rev:3  
2022.02