




**MICRO-TECH**<sup>®</sup>  
ENDOSCOPY

## EFFICIENT RESECTION COLD SNARES

Designed specifically to resect diminutive polyps, the Micro-Tech cold snare is available in two sizes, 10 mm and 15 mm. The thin cutting wire provides for a guillotine-type cut while the loop shape assists in capturing smaller polyps. A stiffer catheter aids in insertion through an extended scope position while providing support during resection. These features combine to provide an effective, efficient and economical solution to resecting diminutive polyps.

OVAL OR DIAMOND LOOP

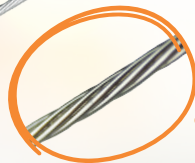

0.23mm WIRE

HANDLE MAY  
BE ROTATED

SMOOTH GLIDE

FIRM CATHETER

NEVER  
COMPROMISE  
ON QUALITY

DRAMATICALLY  
IMPROVE YOUR  
BOTTOM LINE

THROW AWAY  
CONTRACTS  
FOREVER

RELIABLE  
SUPPLY  
PARTNERSHIP

## KEY BENEFITS

### DEDICATED COLD SNARE

Specifically designed for resecting diminutive polyps using a cold snare technique, the oval or diamond shape loop and flexible wire, capture smaller polyps with ease.

### THIN WIRE

The thin 0.23 mm cutting wire provides a guillotine cut that reduces the force needed to complete the resection and may limit tearing or snagging.

### COMFORTABLE ERGONOMIC HANDLES

Accommodating various hand sizes, the ergonomic three-ring handle features a shortened throw and comfortable grip. The stiff catheter provides support for a smooth, effortless cut.

## SPECIFICATIONS

### COLD SNARES

Order Number	Henry Schein Item Number	Description	Loop Diameter (mm)	Sheath Diameter (Fr)	Working Length (cm)	Minimum Channel Size (mm)	Package Units
CS50011	129-0671	Diamond Thin Wire	10	7	230	2.8	10/Box
CS50021	129-0672	Diamond Thin Wire	15	7	230	2.8	10/Box
CS50031	129-0160	Oval Thin Wire	10	7	230	2.8	10/Box
CS50041	129-0161	Oval Thin Wire	15	7	230	2.8	10/Box

## CAN'T FIND IT?

Additional items may be available. Contact us if you can't find what you need.



**EXCEPTIONAL**  
ENDOSCOPIC  
SOLUTIONS.

2855 Boardwalk Drive  
Ann Arbor, MI 48104  
877.552.4027

**MTENDOSCOPY.COM**

CustomerCare@MTEndoscopy.com  
MPW50000 Rev. 3