


## PRECISE PLACEMENT **FULLY COVERED & PARTIALLY COVERED ESOPHAGEAL STENTS**

*Precise stent placement and maintaining position throughout the deployment process are essential.* Micro-Tech Endoscopy has over two decades of experience making high-quality, affordable esophageal stent options.

The deployment system has a pliable and tapered tip to facilitate placement through tight strictures while the pre-loaded delivery catheter and sliding handle mechanism provide smooth delivery of the stent to the targeted site.

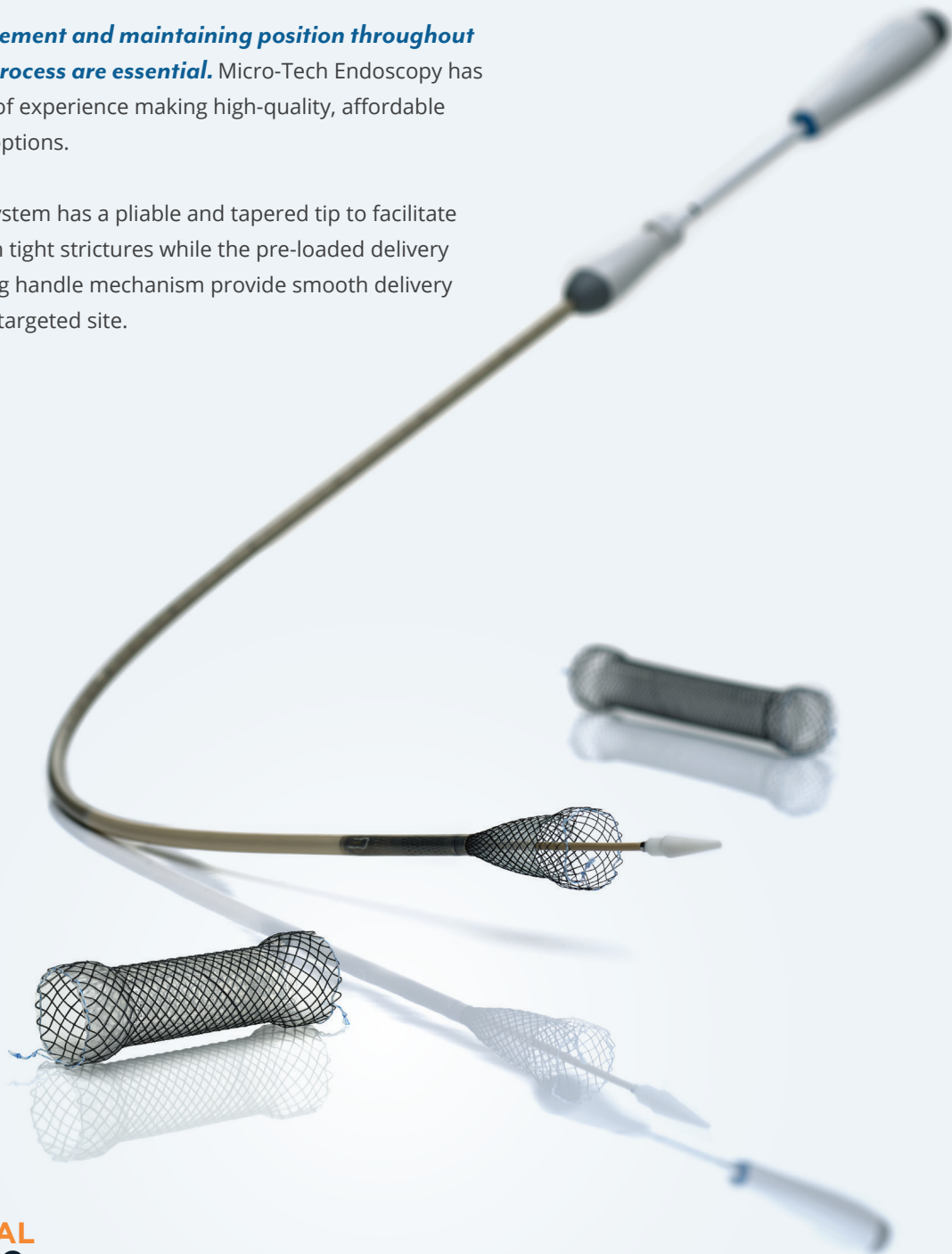

## KEY BENEFITS

### RADIO-OPAQUE MARKERS

Although the entire stent is visible under X-ray fluoroscopy, each stent has several radio-opaque markers to assist with precise stent placement.

### MANUFACTURED FROM ONE WIRE

Our fully and partially covered esophageal stents are braided from a single Nitinol wire, reducing the risk of fracture at the connection points.

### MIGRATION RESISTANCE

The flared ends may reduce migration within the esophageal lumen.


## SPECIFICATIONS

### PARTIALLY COVERED ESOPHAGEAL STENT



Order Number	Description	Diameter of Stent (mm)	Diameter of flares (mm)	Stent Length (cm)	Catheter Outer Diameter (mm)	Working Length (cm)
ET80070	Partially Covered Esophageal Stent	20	25	6	6.7	65
ET80080	Partially Covered Esophageal Stent	20	25	8	6.7	65
ET80090	Partially Covered Esophageal Stent	20	25	10	6.7	65
ET80100	Partially Covered Esophageal Stent	20	25	12	6.7	65
ET80110	Partially Covered Esophageal Stent	20	25	14	6.7	65
ET80120	Partially Covered Esophageal Stent	20	25	15	6.7	65

### FULLY COVERED ESOPHAGEAL STENT

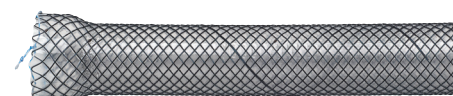

Order Number	Description	Diameter of Stent (mm)	Diameter of flares (mm)	Stent Length (cm)	Catheter Outer Diameter (mm)	Working Length (cm)
ET80190	Fully Covered Esophageal Stent	18	23	6	6.7	65
ET80200	Fully Covered Esophageal Stent	18	23	8	6.7	65
ET80210	Fully Covered Esophageal Stent	18	23	10	6.7	65
ET80220	Fully Covered Esophageal Stent	18	23	12	6.7	65
ET80230	Fully Covered Esophageal Stent	18	23	14	6.7	65
ET80250	Fully Covered Esophageal Stent	20	25	6	6.7	65
ET80260	Fully Covered Esophageal Stent	20	25	8	6.7	65
ET80270	Fully Covered Esophageal Stent	20	25	10	6.7	65
ET80280	Fully Covered Esophageal Stent	20	25	12	6.7	65
ET80290	Fully Covered Esophageal Stent	20	25	14	6.7	65
ET80300	Fully Covered Esophageal Stent	20	25	15	6.7	65
ET80310	Fully Covered Esophageal Stent	22	27	6	6.7	65
ET80320	Fully Covered Esophageal Stent	22	27	8	6.7	65
ET80330	Fully Covered Esophageal Stent	22	27	10	6.7	65
ET80340	Fully Covered Esophageal Stent	22	27	12	6.7	65
ET80350	Fully Covered Esophageal Stent	22	27	14	6.7	65

Note: The inner tube contains a central lumen that accommodates a 0.038 in (0.97 mm) guide wire.

### CAN'T FIND IT?

Additional items may be available. Contact us if you can't find what you need.



2855 Boardwalk Drive  
Ann Arbor, MI 48104  
877.552.4027

[MTENDOSCOPY.COM](http://MTENDOSCOPY.COM)  
customer@mtendoscopy.com

MPW80000 Rev.3