


PRECISE PLACEMENT SEGMENTED ESOPHAGEAL STENTS


Precise stent placement and maintaining position throughout the deployment process are essential.

Micro-Tech Endoscopy has over two decades of experience making high quality, affordable esophageal stent options. Our segmented esophageal stent conforms to bypass esophageal stenoses and allows greater flexibility which may be more comfortable for patients.

The deployment system has a pliable and tapered tip to facilitate placement through tight strictures while the pre-loaded delivery catheter and sliding handle mechanism provide smooth delivery of the stent to the targeted site. Our line of esophageal stents combine flexibility and consistent radial force to help with stricture resolution.

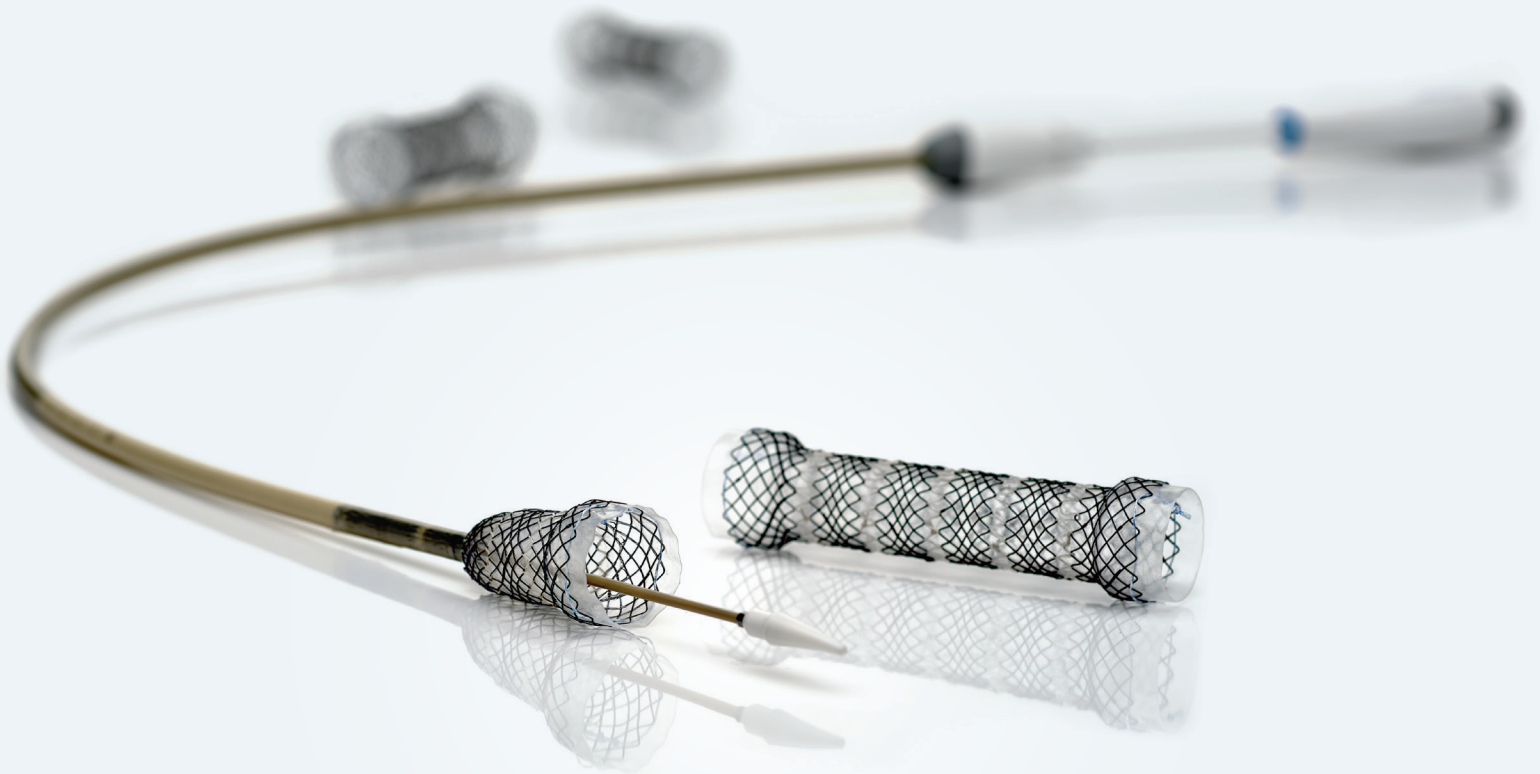

KEY BENEFITS

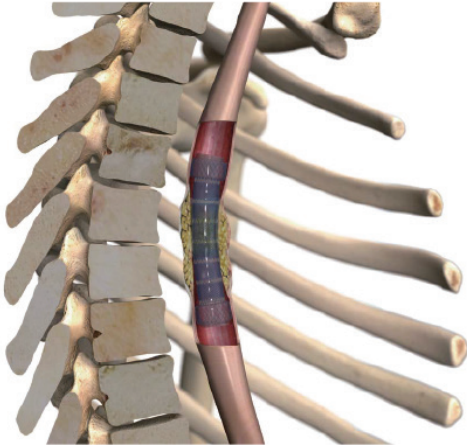

RADIO-OPAQUE MARKER

Although the entire stent is visible under X-ray fluoroscopy, each stent has several radio-opaque markers to assist with precise stent placement.

MIGRATION RESISTANCE

The flared ends along with the unique segmented construction may reduce migration within the esophageal lumen.

INDEPENDENT MOVEMENT

The segmented modular stent allows independent movement between each section which maximizes conformance to curved or tortuous anatomy.

REPOSITIONING SUTURES

Sutures are located on the distal and proximal ends of the stent to allow repositioning or removal immediately post stent placement.

SPECIFICATIONS



SEGMENTED ESOPHAGEAL STENT

Order Number	Description	Diameter of Stent (mm)	Diameter of flares (mm)	Stent Length (cm)	Catheter Outer Diameter (mm)	Working Length (cm)
ET80510	Segmented Esophageal Stent	18	24	6	6.7	65
ET80520	Segmented Esophageal Stent	18	24	8	6.7	65
ET80530	Segmented Esophageal Stent	18	24	10	6.7	65
ET80540	Segmented Esophageal Stent	18	24	12	6.7	65
ET80550	Segmented Esophageal Stent	18	24	14	6.7	65
ET80560	Segmented Esophageal Stent	20	26	6	6.7	65
ET80570	Segmented Esophageal Stent	20	26	8	6.7	65
ET80580	Segmented Esophageal Stent	20	26	10	6.7	65
ET80590	Segmented Esophageal Stent	20	26	12	6.7	65
ET80600	Segmented Esophageal Stent	20	26	14	6.7	65
ET80610	Segmented Esophageal Stent	22	28	6	6.7	65
ET80620	Segmented Esophageal Stent	22	28	8	6.7	65
ET80630	Segmented Esophageal Stent	22	28	10	6.7	65
ET80640	Segmented Esophageal Stent	22	28	12	6.7	65
ET80650	Segmented Esophageal Stent	22	28	14	6.7	65

Note: The inner tube contains a central lumen that accommodates a 0.038 in (0.97 mm) guide wire.

CAN'T FIND IT?

Additional items may be available. Contact us if you can't find what you need.



2855 Boardwalk Drive
Ann Arbor, MI 48104
877.552.4027

MTENDOSCOPY.COM
customer@mtendoscopy.com

MPW85000 Rev.3