


DEPENDABLE & PRECISE ENDOSCOPIC ULTRASOUND NEEDLES

Whether requiring an EUS solution for diagnostic and therapeutic intervention or for histological sampling, Micro-Tech Endoscopy provides high quality, affordable clinical options to help with consistent sampling and accurate diagnosis from EUS tissue acquisition.

The Trident™ FNB Needle with its unique tri-tip core needle design achieves maximized yield and enables consistent histological tissue acquisition.

The Areus™ FNA Needle features an elongated dual edged bevel tip to aid in precise penetration into the target area delivering a consistent yield for cytology.

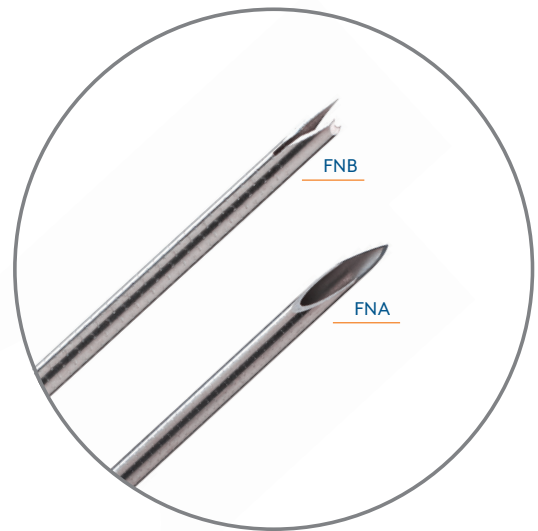

KEY BENEFITS



STYLET ETCHING

Engraving on stylet provides tactile feedback when the stylet is close to exiting the needle shaft. This indicator may help prevent the stylet from releasing unexpectedly resulting in contamination or spray.



RAPID SHEATH AND NEEDLE ADJUSTORS

The push button adjustors on Micro-Tech EUS needles offer quick and secure changes in needle depth and sheath length. The ergonomic handle design with smooth slide action provides precise one-handed control.



STYLET MANAGEMENT CAP

Integrated stylet cap that holds looped stylet securely, so it is easier to manage.

NEEDLE ECHOGENICITY

It is always important to maintain tip visibility for accurate guidance within the target area. The needle surface features a laser etched V pattern, that provides good visibility under ultrasound.

COBALT-CHROMIUM (COCR)

The Trident and Areus EUS needles are made of CoCr alloy which maintains its shape and sharpness while providing durability, flexibility and reliable target accessibility.

SPECIFICATIONS

TRIDENT ENDOSCOPIC ULTRASOUND FINE NEEDLE BIOPSY STERILE

Order Number	Description	Needle Size	Sheath Diameter (mm)	Minimum Channel Size (mm)	Package Units
ES34611	Trident FNB Needle	25ga	1.8	2.8	Each
ES34591	Trident FNB Needle	22ga	1.8	2.8	Each

AREUS ENDOSCOPIC ULTRASOUND ASPIRATION NEEDLE STERILE

Order Number	Description	Needle Size	Sheath Diameter (mm)	Minimum Channel Size (mm)	Package Units
ES34621	Areus FNA Needle	25ga	1.8	2.8	Each
ES34601	Areus FNA Needle	22ga	1.8	2.8	Each

- Working Length: 137.5cm to 141.5cm, adjustable
- Needle Length: 0 to 8cm, adjustable
- Includes Vacuum Lock 20cc syringe and one-way stopcock

CAN'T FIND IT?

Additional items may be available. Contact us if you can't find what you need.



2855 Boardwalk Drive
Ann Arbor, MI 48104
877.552.4027

MTEndoscopy.com
customer@mtendoscopy.com

EUS_MPW34000 Rev.4