


PRECISE CUT HOT SNARES

Our line of hot and cold snares features effective cutting wires and user-friendly handles at cost-effective prices.

With the innovative power connector, the snare can be operated and rotated with ease.

Our hot snares are available in a compliant braided wire or a firmer wire for effective tissue acquisition and resection.

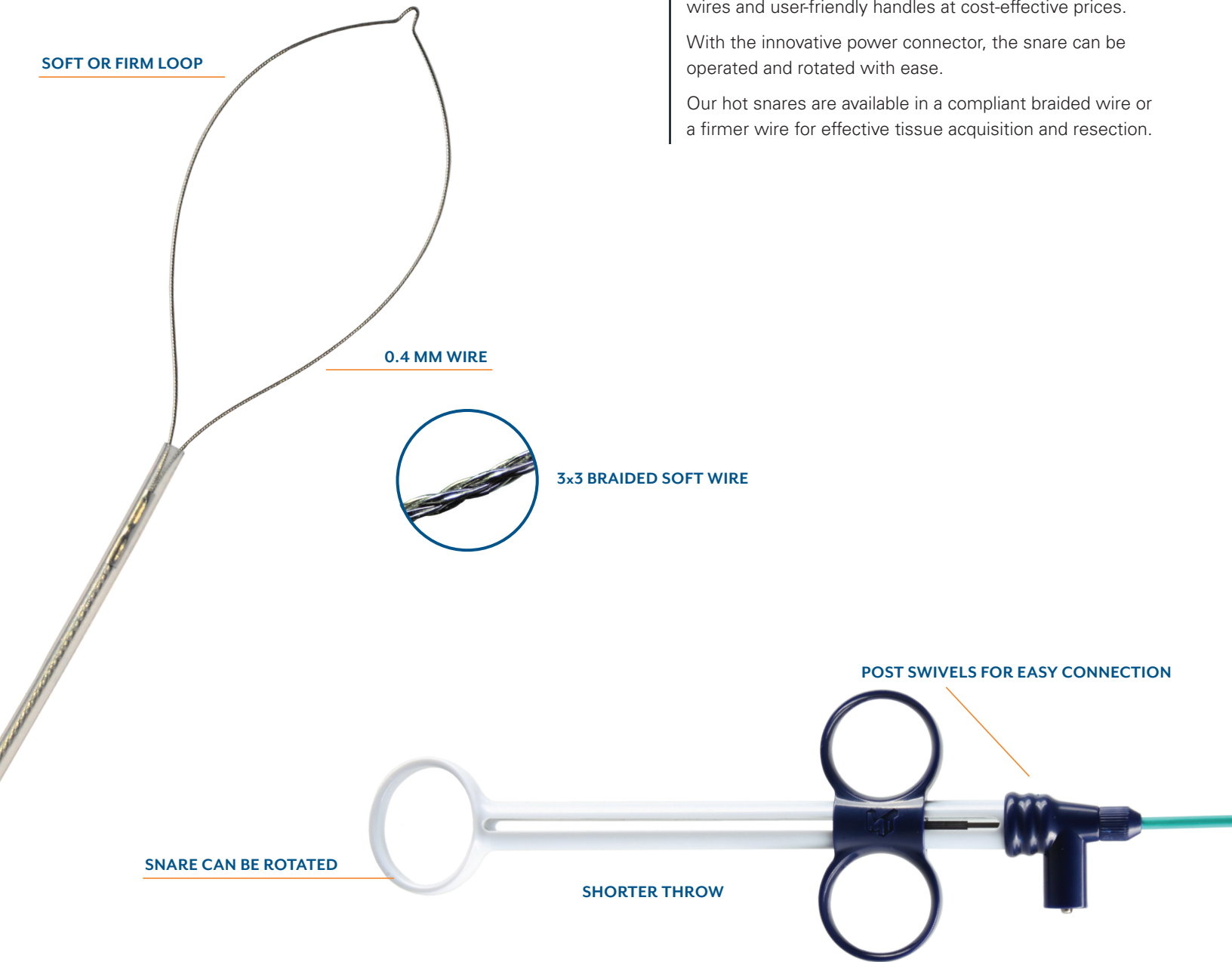

KEY BENEFITS

SHORTER THROW HANDLE

The handle throw is adjusted according to the snare size. Small and Medium snares have a shorter throw than Large or Extra Large snares. A shorter throw improves comfort for users with smaller hands.

BRAIDED WIRE

The braided wire combines the tissue acquisition features of a firm snare with the flexibility and smooth deployment of a soft snare. Suitable for daily use in facilities looking to minimize inventory.

INNOVATIVE POWER CONNECTION

The handle features a dedicated power connection, separate from the finger ring slider. This feature permits unrestricted movement and rotation of the snare handle and may reduce procedural delays.

HOT SNARES BRAIDED STERILE Ⓢ

Order Number	Henry Schein Item Number	Description	Loop Diameter (mm)	Sheath Diameter (Fr)	Working Length (cm)	Minimum Channel Size (mm)	Package Units
PS51011	129-0186	Mini Braided Oval	10	7	230	2.8	10/Box
PS51021	129-0673	Small Braided Oval	15	7	230	2.8	10/Box
PS51031	129-0674	Medium Braided Oval	20	7	230	2.8	10/Box
PS51041	129-0675	Large Braided Oval	24	7	230	2.8	10/Box
PS51051	129-0187	Extra Large Braided Oval	30	7	230	2.8	10/Box

HOT SNARES FIRM STERILE Ⓢ

Order Number	Henry Schein Item Number	Description	Loop Diameter (mm)	Sheath Diameter (Fr)	Working Length (cm)	Minimum Channel Size (mm)	Package Units
PS51071	129-0188	Mini Oval Firm	10	7	230	2.8	10/Box
PS51081	129-0189	Small Oval Firm	15	7	230	2.8	10/Box
PS51091	129-0190	Medium Oval Firm	20	7	230	2.8	10/Box
PS51101	129-0191	Large Oval Firm	24	7	230	2.8	10/Box
PS51111	129-0192	Extra Large Oval Firm	30	7	230	2.8	10/Box
PS51121	129-0193	Jumbo Oval Firm	36	7	230	2.8	10/Box

HOT SNARES HEX FIRM STERILE Ⓢ

Order Number	Henry Schein Item Number	Description	Loop Diameter (mm)	Sheath Diameter (Fr)	Working Length (cm)	Minimum Channel Size (mm)	Package Units
PS52031	129-0194	Small Hex Firm	15	7	230	2.8	10/Box
PS52041	129-0195	Large Hex Firm	25	7	230	2.8	10/Box

CAN'T FIND IT?

Additional items may be available. Contact us if you can't find what you need.



2855 Boardwalk Drive
Ann Arbor, MI 48104
877.552.4027

MTEndoscopy.com
customer@mtendoscopy.com

MPW51000 Rev. 4