


PRECISION CUT LESIONHUNTER COLD SNARE



LesionHunter™ Cold Snare features a rotatable thin Nitinol Wire which maintains the oval shape cut after cut. With this cold snare, you get a clean and precise cut.

- Dedicated Cold Snare – Specifically designed for resecting polyps using a cold snare technique.
- Thin Nitinol Wire – This wire material provides durability in the shape and allows for a 0.18 mm thin wire for an even cleaner cut.
- Coiled Sheath – The coiled metal catheter creates greater stability for the guillotine-like cutting action of the wire.
- Comfortable Ergonomic Handles – Like other Micro-Tech snares, this snare has a three-ring comfortable grip with shortened throw for an effortless cut.
- Reduced Fly-Away Tissue - The thin wire may reduce the risk of generating a “fly-away” specimen; thereby making it easier to find the tissue at the end of the resection.

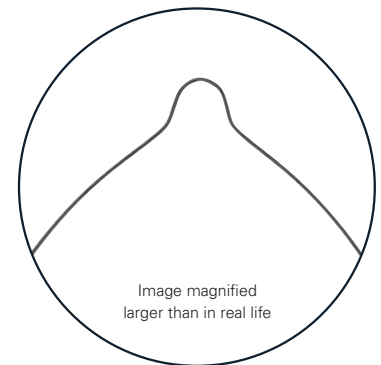

SPECIFICATIONS

LESIONHUNTER STERILE

Order Number	Description	Length (cm)	Outer Diameter (mm)	Package Units
CS50051	LesionHunter Cold Snare, Oval 10 mm wide, Nitinol Wire	230	2.3	10/Box
CS50061	LesionHunter Cold Snare, Oval 15 mm wide, Nitinol Wire	230	2.3	10/Box
CS50071	LesionHunter Cold Snare, Oval 20 mm wide, Nitinol Wire	230	2.3	10/Box

CAN'T FIND IT?

Additional items may be available. Contact us if you can't find what you need.