


EFFECTIVE LIGATION MULTIPLE BAND LIGATORS



Endoscopic banding has been demonstrated as an effective and simple treatment for internal hemorrhoids and varices. Hemorrhoid banding is readily performed in an ambulatory setting and can result in less pain and have a shorter recovery period than a surgical hemorrhoidectomy.

The Multi-Band Ligator Sets from Micro-Tech offer cost-effective solutions for the treatment of hemorrhoids and varices during routine endoscopy procedures.

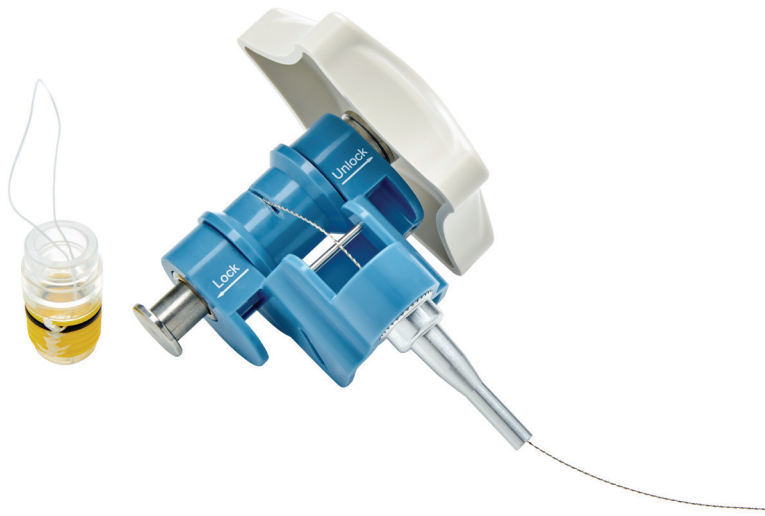

KEY BENEFITS

EXCELLENT VISIBILITY

A clear ligating barrel with rear mounted bands optimizes the endoscopic view while maintaining suction.

BAND DESIGN

The square band design provides consistent shape and reliable retention. The hemorrhoid set comes with four pre-loaded latex-free bands for routine colonoscopies. The upper ligator comes pre-loaded with six latex-free bands and can be used in upper and lower endoscopic procedures.

USER-FRIENDLY SYSTEM

The handle unit comes with the pull wire pre-loaded to simplify setup. The alert band notifies the operator when only one band remains.

SPECIFICATIONS

MULTI-BAND LIGATOR STERILE

Order Number	Description	Minimum Working Channel (mm)	Working Length (cm)	Package Units
HL32161	Endoscopic Hemorrhoid Multi-Band Ligator Set, 4 Latex Free Bands, compatible with 9.4-13 mm scope diameter	2.8	220	1/Box
ML32651	Endoscopic Multi-Band Ligator Set, 6 Latex Free Bands, compatible with 9.4-13 mm scope diameter	2.8	145	1/Box

CAN'T FIND IT?

Additional items may be available. Contact us if you can't find what you need.