


## IDEAL BITE BIOPSY FORCEPS

With our wide selection of cost-effective disposable biopsy forceps, you are perfectly equipped for every procedure.

Available jaw diameters of 1.8 mm, 2.3 mm and 3.0 mm, and are offered in Oval, Serrated, Alligator jaws, with or without a spike.

- Available lengths of 160 cm, 180 cm and 230 cm
- All forceps feature entry and exit indicator markings

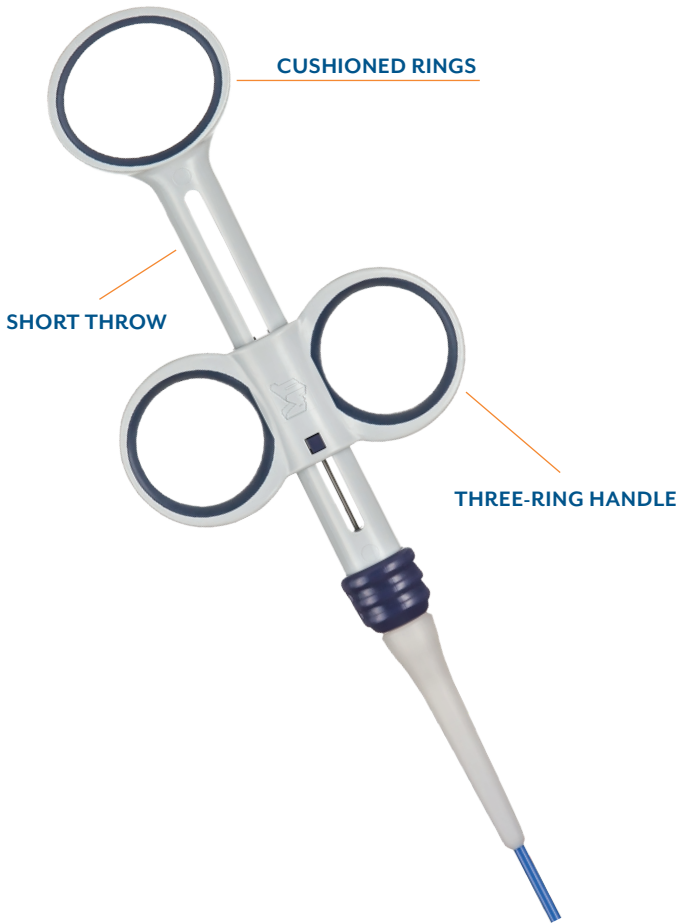

### KEY BENEFITS

#### ERGONOMIC HANDLE

Three-ring comfort with cushioned lining and short throw provides improved grip and accommodates smaller hands.

#### DIRECT DRIVE SYSTEM

The single drive wire and hinge linkage deliver smooth activation around curves and tortuous anatomy. The system delivers more closure force to the jaws which provides for improved cutting ability and assures complete jaw closure.

#### SCOPE FRIENDLY DESIGN

The design of the Micro-Tech forceps family is based upon the robust and proven design of original reusable biopsy instruments. These designs provide optimal bite, while minimizing the risk of scope damage.

## BIOPSY FORCEP FOR 3.2 MM SCOPE CHANNELS STERILE ⌘

Order Number	Description	Jaw Style	Sheath Color	Working Length (cm)	Minimum Channel Size (mm)	Package Units
BF40061	Jumbo Biopsy Forceps	Oval	Blue	230	3.2	10/Box
BF40071	Jumbo Biopsy Forceps with Spike	Oval	Blue	230	3.2	10/Box
BF40171	Jumbo Biopsy Forceps	Alligator	Blue	230	3.2	10/Box
BF40181	Jumbo Biopsy Forceps with Spike	Alligator	Blue	230	3.2	10/Box

## BIOPSY FORCEP FOR 2.8 MM SCOPE CHANNELS STERILE ⌘

Order Number	Description	Jaw Style	Sheath Color	Working Length (cm)	Minimum Channel Size (mm)	Package Units
BF40081	Standard Biopsy Forceps	Oval	Orange	180	2.8	10/Box
BF40111	Standard Biopsy Forceps with Spike	Oval	Orange	180	2.8	10/Box
BF40091	Standard Biopsy Forceps	Oval	Blue	230	2.8	10/Box
BF40121	Standard Biopsy Forceps with Spike	Oval	Blue	230	2.8	10/Box
BF40096	Standard Biopsy Forceps	Oval	Blue	230	2.8	50/Box
BF40126	Standard Biopsy Forceps with Spike	Oval	Blue	230	2.8	50/Box
BF40191	Standard Biopsy Forceps	Alligator	Orange	180	2.8	10/Box
BF40221	Standard Biopsy Forceps with Spike	Alligator	Orange	180	2.8	10/Box
BF40561	Standard Biopsy Forceps	Alligator	Blue	230	2.8	10/Box
BF40541	Standard Biopsy Forceps with Spike	Alligator	Blue	230	2.8	10/Box
BF40566	Standard Biopsy Forceps	Alligator	Blue	230	2.8	50/Box
BF40546	Standard Biopsy Forceps with Spike	Alligator	Blue	230	2.8	50/Box
BF40431	Standard Biopsy Forceps	Serrated	Orange	180	2.8	10/Box
BF40411	Standard Biopsy Forceps with Spike	Serrated	Orange	180	2.8	10/Box
BF40441	Standard Biopsy Forceps	Serrated	Blue	230	2.8	10/Box
BF40421	Standard Biopsy Forceps with Spike	Serrated	Blue	230	2.8	10/Box
BF40446	Standard Biopsy Forceps	Serrated	Blue	230	2.8	50/Box
BF40426	Standard Biopsy Forceps with Spike	Serrated	Blue	230	2.8	50/Box

## BIOPSY FORCEP HOT STERILE ⌘

Order Number	Description	Jaw Style	Sheath Color	Working Length (cm)	Minimum Channel Size (mm)	Package Units
HF48521	Hot Biopsy Forceps Lower	Oval	Red	230	2.8	10/Box
HF48541	Hot Biopsy Forceps Lower	Serrated	Red	230	2.8	10/Box

## ACTIVE CORD

Order Number	Description	Package Units
QB51011	3.0 mm Active Cord, Reusable, Bovie 8.0 mm Male Fitting	1/Box

## CAN'T FIND IT?

Additional items may be available. Contact us if you can't find what you need.



2855 Boardwalk Drive  
Ann Arbor, MI 48104  
877.552.4027

[MTEndoscopy.com](http://MTEndoscopy.com)  
[customer@mtendoscopy.com](mailto:customer@mtendoscopy.com)

MPW40000 Rev. 6