


## PRECISE CANNULATION SPHINCTEROTOME

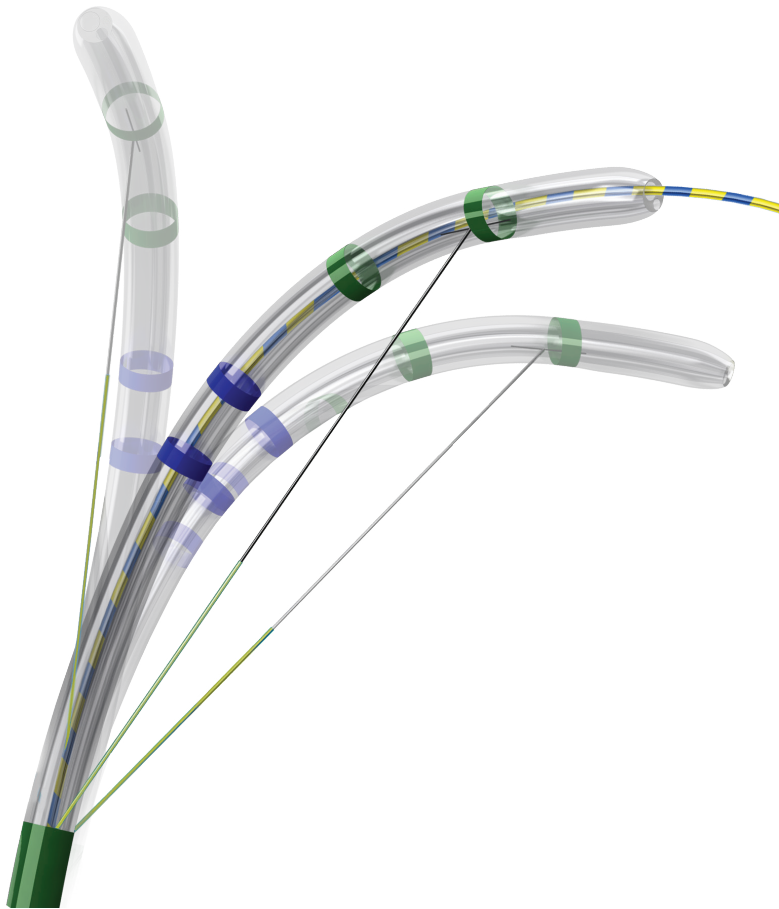

The Micro-Tech sphincterotome features a patent pending\* opening mechanism to facilitate both short wire and long wire procedures with one device. The rotatable design and locking mechanism provide maximum control during cannulation and sphincterotomy.

### KEY BENEFITS

#### FLEXIBILITY

Compatible with wires up to .035", the rotatable tip can be actuated to ensure proper cannulation and sphincterotomy positioning.

#### CONTROL

The c-channel catheter allows for a splitting movement to separate wire from catheter for physician guided cannulation.

#### ATRAUMATIC TIP


The tapered atraumatic tip helps to ease insertion and tracking.

\*US patent pending

# MAXIMUM FLEXIBILITY

## BILIARY STONE RETRIEVAL

Micro-Tech's stone extraction balloon is comprised of a high-grade latex for reliability during endoscopic bile stone extraction. The 3 stage balloons come in above and below injection options to meet a wide range of clinical needs. The catheter features a 3-lumen design for insufflation, guidewire access, and contrast injection.

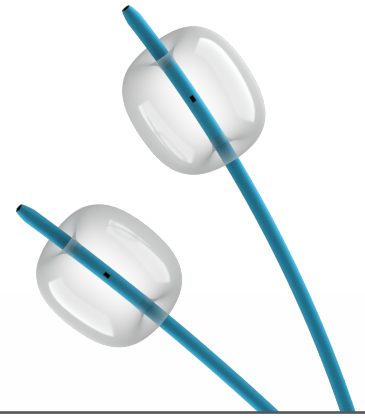

### KEY BENEFITS

#### PROCEDURAL FLEXIBILITY

3-stage balloons are sized in 3 mm increments between 9 and 21 mm.

#### INSERTION AND TRACKING

The tapered atraumatic tip helps to ease insertion and improve tracking along the wire into biliary anatomy.

#### COMPATIBILITY

The balloon catheter is compatible with guidewires up to .035".

#### SPHINCTEROTOMES STERILE STERILE EO

Order Number	Tip Length (mm)	Cutting Wire Insulation	Cut Wire	Preloaded Wire Length (cm)	Working Length (cm)	Package Units
SS75041	5	Insulated	20 mm	260	200	1/Box
SS75061	5	Insulated	20 mm	450	200	1/Box
SS75101	5	Insulated	25 mm	260	200	1/Box
SS75161	5	Insulated	30 mm	260	200	1/Box

Preloaded with 0.035" Hydro-Slide Straight Tip Guidewire.

#### RETRIEVAL BALLOONS STERILE STERILE EO

Order Number	Catheter Diameter (Fr)	Balloon Stages	Injection Port Location	Working Length (cm)	Package Units
SB75921	7	9/12/15 mm	Above (distal to balloon)	200	1/Box
SB75931	7	12/15/18 mm	Above (distal to balloon)	200	1/Box
SB75941	7	15/18/21 mm	Above (distal to balloon)	200	1/Box
SB76131	7	9/12/15 mm	Below (proximal to balloon)	200	1/Box
SB76141	7	12/15/18 mm	Below (proximal to balloon)	200	1/Box
SB76151	7	15/18/21 mm	Below (proximal to balloon)	200	1/Box

#### CAN'T FIND IT?

Additional items may be available. Contact us if you can't find what you need.



2855 Boardwalk Drive  
Ann Arbor, MI 48104  
877.552.4027

[MTEndoscopy.com](http://MTEndoscopy.com)  
customer@mtendoscopy.com

MPW75000 Rev. 1