


# LOCKADO HEMOSTASIS CLIP



## LOCK INTO TISSUE

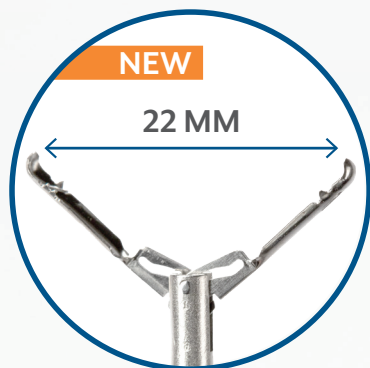

**Lockado<sup>®</sup> allows unrestricted opening and closing so it can be repositioned until the desired position is achieved.**

Hemostasis clips need to provide accuracy, reliability and precise functionality. With two sets of opposing atraumatic teeth in the arms of the clip, Lockado delivers a reliable and secure grip for consistent performance. Driven by Micro-Tech Endoscopy's proven clip delivery system, the Lockado clip provides exceptional reliability and tissue management.


KEY FEATURES

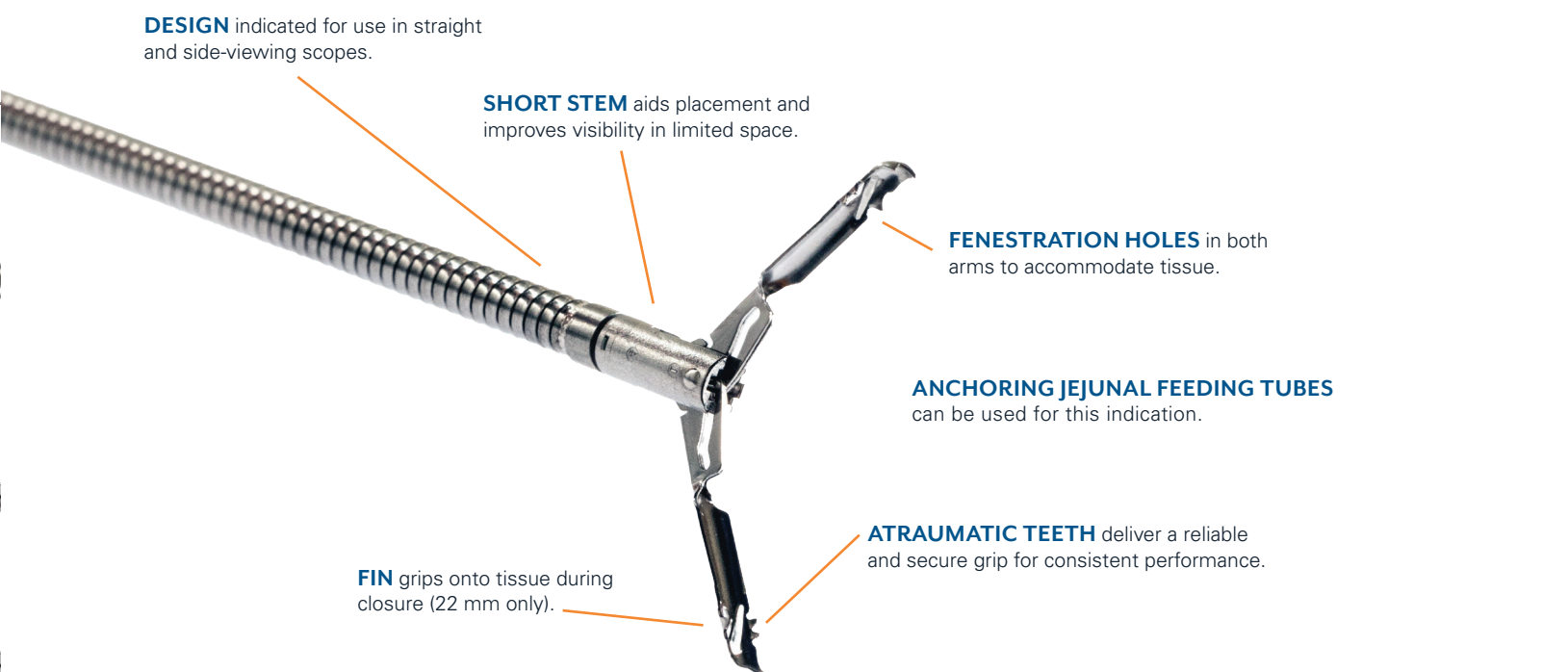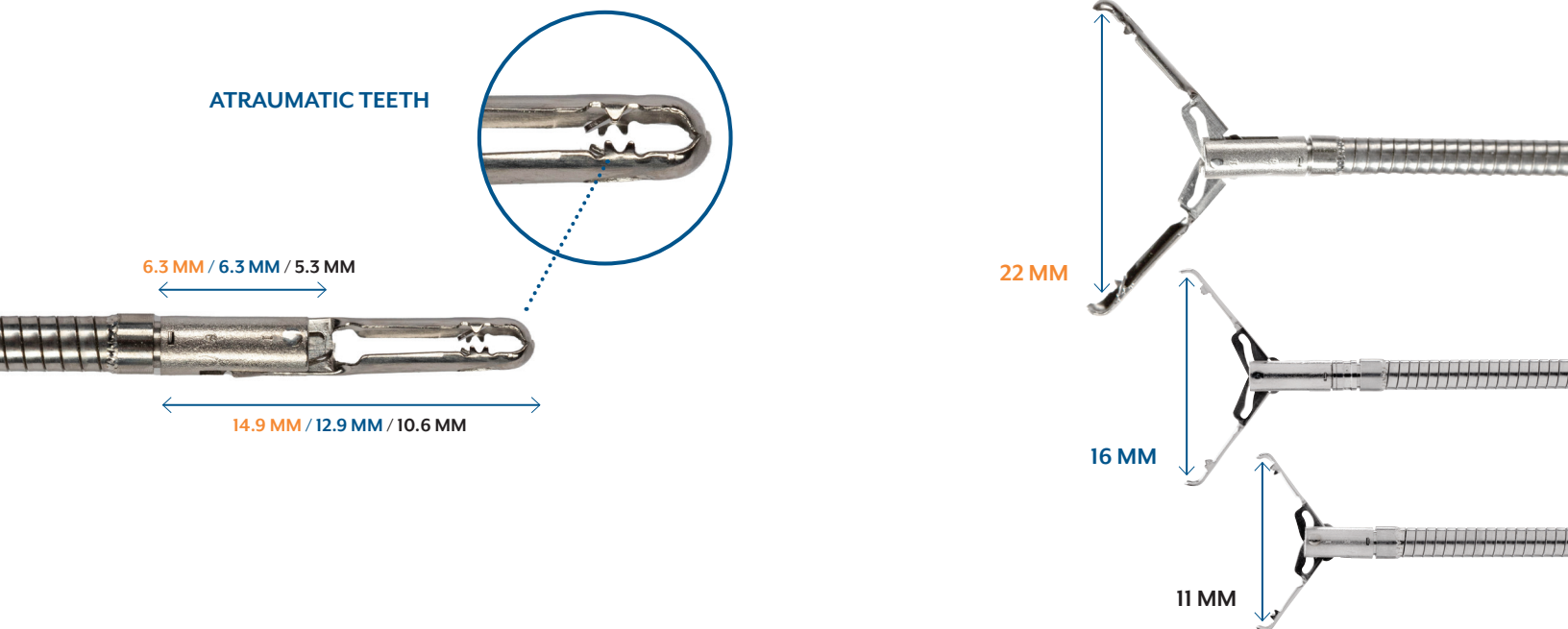


“I prefer the Micro-Tech clips in situations such as post-ERCP sphincterotomy/ampullectomy bleed, perforation closure, EMR/ESD, mucosectomy in esophagus, small bowel or colon because they grip tissue easily and the exceptionally short stems are great in tight spaces.”

- DR. FENG LI, MD, ADVANCED ENDOSCOPIST

**EXCELLENT RETENTION STRENGTH**  
When comparing Lockado to other competitive hemostasis clips, testing showed that all sizes of Lockado have greater retention strength on average than Boston Scientific Resolution® 11 mm and Resolution Ultra 360™.

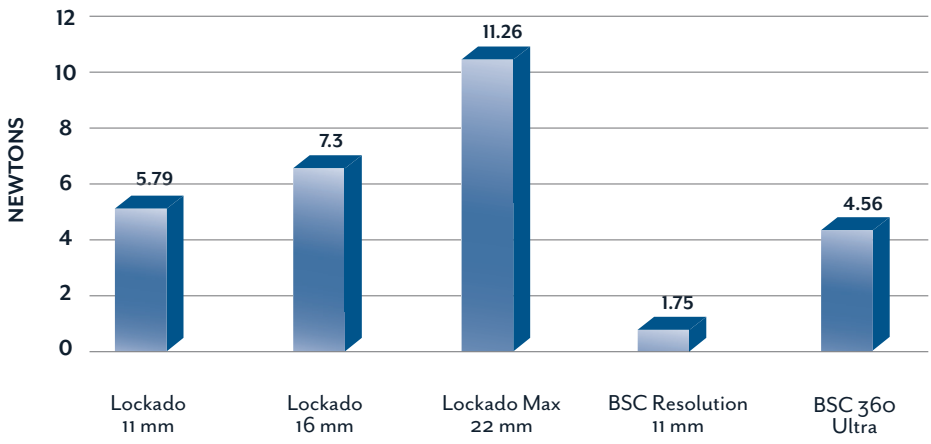

\*Measurements represent average force required to pull clips off simulated tissue.

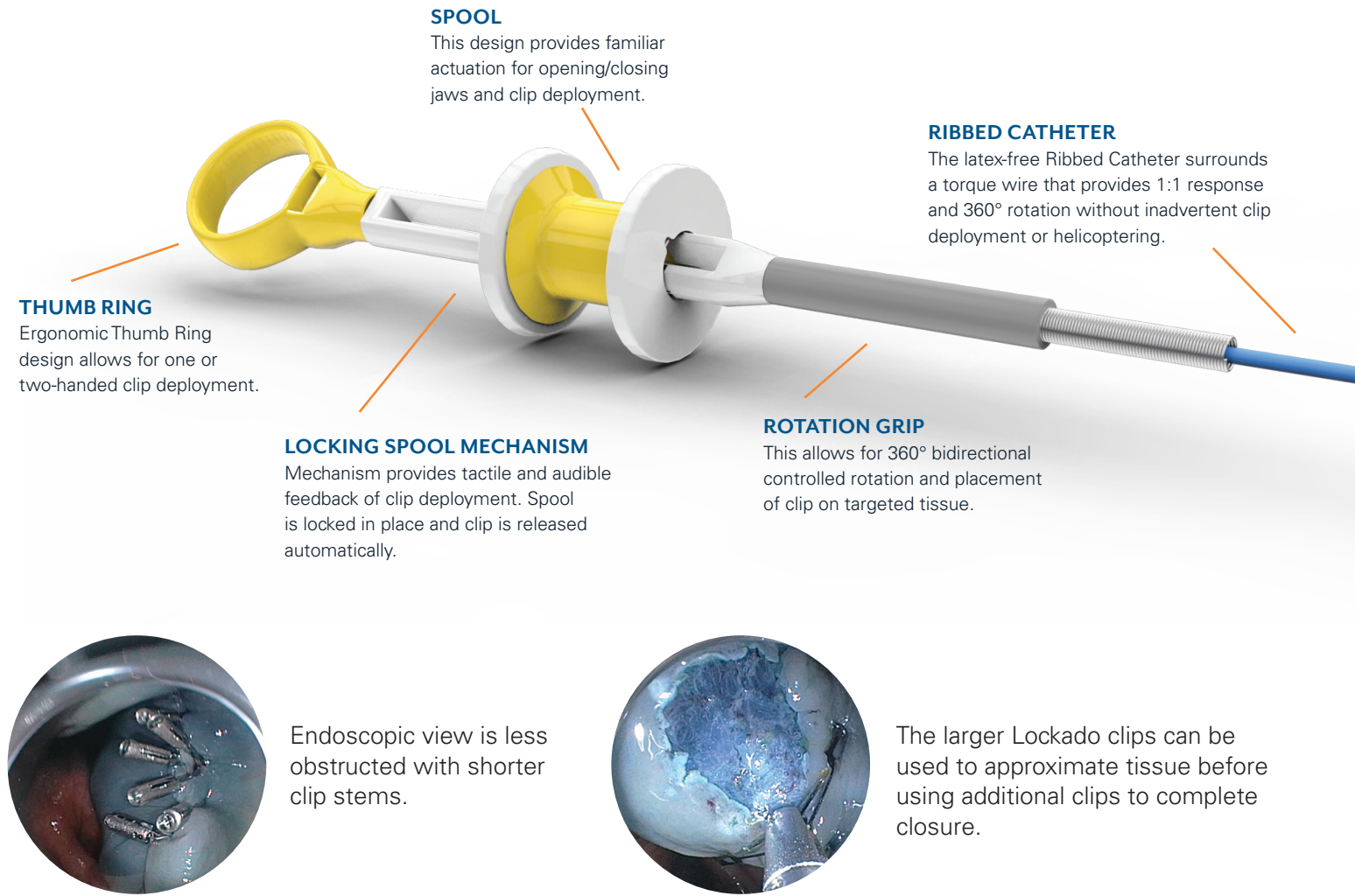

## SPECIFICATIONS

### LOCKADO HEMOSTASIS CLIP STERILE



Order Number	Opening Width (mm)	Sheath Diameter (mm)	Working Length (cm)	Minimum Channel Size (mm)	Packing Units
LO30841	11	2.6	165	2.8	10/Box
LO30861	11	2.6	235	2.8	10/Box
LO30865	11	2.6	235	2.8	2/Box

### LOCKADO PLUS HEMOSTASIS CLIP STERILE



Order Number	Opening Width (mm)	Sheath Diameter (mm)	Working Length (cm)	Minimum Channel Size (mm)	Packing Units
LO30871	16	2.6	165	2.8	10/Box
LO30891	16	2.6	235	2.8	10/Box
LO30895	16	2.6	235	2.8	2/Box

### LOCKADO MAX HEMOSTASIS CLIP\* STERILE



Order Number	Opening Width (mm)	Sheath Diameter (mm)	Working Length (cm)	Minimum Channel Size (mm)	Packing Units
LO35021	22	2.6	165	2.8	10/Box
LO35041	22	2.6	235	2.8	10/Box
LO35045	22	2.6	235	2.8	2/Box

\*Jaws may not open to maximum reach until placed in contact with tissue.

## INDICATIONS FOR USE

The Lockado Hemostasis Clip is indicated for Endoscopic clip placement within the Gastrointestinal tract of adult patients only via a straight or side-viewing flexible endoscope for the purpose of:

- (1) Endoscopic marking;
- (2) Hemostasis for
  - (a) Mucosal / sub-mucosal defects < 3 cm,
  - (b) Bleeding ulcers,
  - (c) Polyps < 1.5 cm in diameter,
  - (d) Diverticula in the colon,
  - (e) Arteries < 2 mm
  - (f) Prophylactic clipping to reduce the risk of delayed bleeding post lesion resection;
- (3) As a supplementary method, closure of GI tract luminal perforations <20 mm that can be treated conservatively.

## MR CONDITIONAL

Non-clinical testing has demonstrated the Lockado Hemostasis Clip is MR Conditional according to ASTM F2503.

<sup>1</sup> Bench testing was performed (n=120) per device. Test report: MTN-SD-TR-001.

<sup>2</sup> Bench testing was performed (n=5, n=25) per device. Test report: MTN-0035150 A1 Comparison Test Report.

All trademarks are the property of their respective owners.

### CAN'T FIND IT?

Additional items may be available. Contact us if you can't find what you need.



2855 Boardwalk Drive  
Ann Arbor, MI 48104  
877.552.4027

[MTEndoscopy.com](http://MTEndoscopy.com)  
[customercare@mtendoscopy.com](mailto:customercare@mtendoscopy.com)

MPW30800B Rev. 2