


SIMPLE & EFFECTIVE FOR  
LARGE STONE MANAGEMENT

# EHL LITHOTRON EL27 AND ACCESSORIES

The Walz Lithotron EL27 Compact offers a secure, easy, and convenient way to administer bipolar ElectroHydraulic Lithotripsy (EHL).

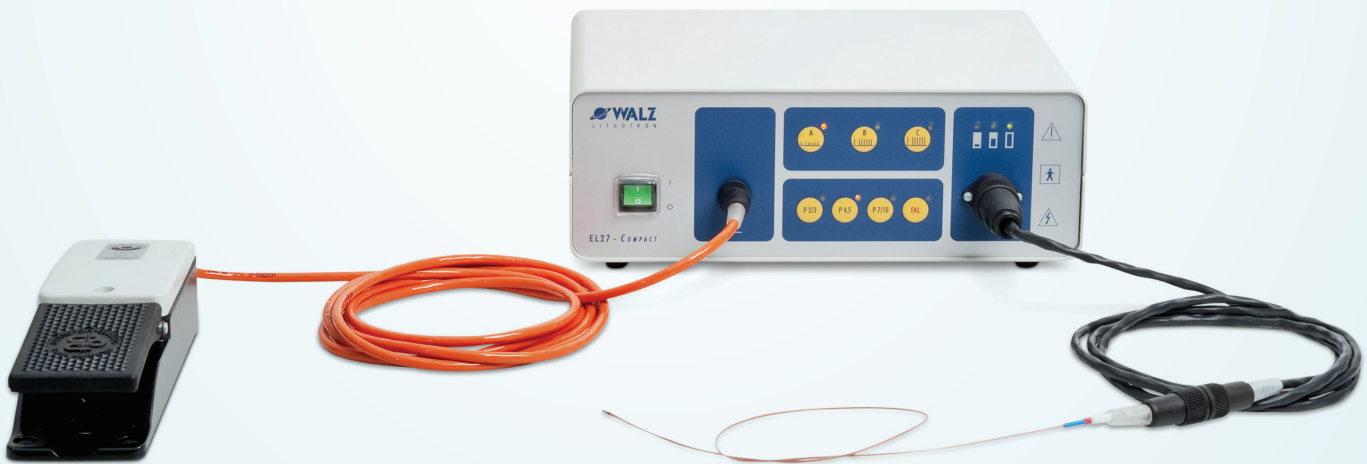

## KEY BENEFITS

### CONTROL

Offering 3 different intensity options (up to 900 mJ), allowing users the ability to adjust energy levels based on stone densities.

### SECURITY

Reduced tensile phases of concentrated energy enable accurate stone fragmentation while minimizing harm to soft tissues.


### EFFICIENCY

The centered micro-shockwaves delivered by the probe facilitates targeting with fast pulses and sharply rising amplitudes.

## EHL PROBES

- Highly flexible with enough tensile strength to maintain position.
- Precision machined probe tip for reproduceable shockwaves.
- Rounded probe tip edges for safe and smooth passage through endoscope.
- Centered micro spark plasmas preserve distances to the tissue and facilitates targeting.
- Multiple probe size options aiding in increased surface area contact with stone for efficient disintegration of concrements.

### Electrohydraulic Lithotripsy (EHL)

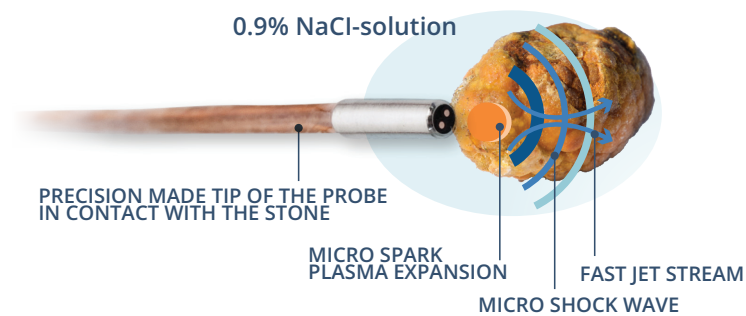

Order Number	Probe Diameter (Fr)	Probe Diameter (mm)	Working Length (cm)
P243000FST	2.4	0.8	300
P303000FST	3.0	1.0	300
P453000FST	4.5	1.5	300

\*\*Note: 3.0 Fr and 4.5 Fr EHL probes to be used exclusively with 3.9mm eyeMax scope with 2.0mm working channel.

#### CAN'T FIND IT?

Additional items may be available. Contact us if you can't find what you need.

